

CLINICALLY PROVEN AND
BREATHABLE HIP PROTECTION

SAFEHIP[®] AIR X



ABOUT SAFEHIP AIR X ↘

SAFEHIP[®] AirX is a 100% textile hip protector with the protector shield based on AirX spacer fabric technology, which is a breathable, skin friendly and soft textile material.

The SAFEHIP[®] AirX hip protectors are based on our specially horseshoe-shaped textile shield. The shield technology of SAFEHIP[®] uses an energy shunting technique to reduce the fall impact by dispersing the energy from the fall away from the hip area. SAFEHIP[®] AirX is an evidence-based, clinical and very comfortable hip fracture prevention solution intended to reduce the incidence of hip fractures.

FEATURES

- Soft and comfortable to wear, seamlessly knitted pants
- Easy to take on and off
- 100% textile all the way through
- 100% breathable
- Minimises moisture build-up and feels comfortable when in contact with the skin
- Considerably more comfortable than conventional hip protectors with rigid non-textile shields
- Designed to disperse energy from a fall away from the hip joint and to resist multiple falls
- The shields adapt to the body perfectly
- Seamless knitting with the purpose of avoiding skin discomfort for highest user comfort
- Suitable for incontinent users



HIP FRACTURE PREVENTION



SEAMLESS KNITTING



CLINICALLY PROVEN



FACT-BASED SIZING

EUROPEAN-BASED PRODUCTION
FOR OPTIMAL QUALITY,
SUPPLY CHAIN EFFICIENCY
AND FLEXIBILITY



STANDARD
100
6234-7001
DTI

PRODUCT RANGE

SAFEHIP® AIR X



SAFEHIP® AIRX UNISEX SEWN-IN OR REMOVABLE 16 MM SHIELDS, WHITE

Art. no. 336550 (sewn-in)
Art. no. 336551 (removable)

SAFEHIP® AIRX FEMALE SEWN-IN 16 MM SHIELDS, WHITE

Art. no. 3366

SAFEHIP® AIRX MALE SEWN-IN 16 MM SHIELDS, WHITE

Art. no. 3367

SAFEHIP® AIRX OPEN SEWN-IN 16 MM SHIELDS, WHITE

Art. no. 3368

PRODUCT DESCRIPTION

The designs of SAFEHIP® AirX hip protector pants ensure optimum protection. The hip protectors come in four models – a female, male, unisex and an open model in white. The pants are flexible and comfortable to wear and may be used with incontinence pads as a fixation pant.

The male, female and unisex versions have been designed to perform like regular underwear, with a fly opening on the male model and lace-look detail on the female. The Unisex brief is suitable for both.

The open model offers a great degree of protection for people suffering from heavy or double incontinence. Designed to be worn underneath an adult diaper, the large crutch opening makes it quick and easy to change diapers.

CLINICAL EVIDENCE

SAFEHIP® is the most thoroughly tested hip protector in the world. It has been tested on more than 7,000 users in clinical trials in Scandinavia, Holland, Germany, Japan and Australia, and the results are unequivocal: SAFEHIP® has a significant effect in reducing the number of hip fractures.

BETTER USER COMFORT AND COMPLIANCE WITH THE AIR X® TEXTILE SHIELD


The SAFEHIP® AirX hip protectors feature a textile shield made of AirX which is a 3D mesh fabric, ideally suited to absorb impacts and cushion the body from falls. The material has the ability to efficiently transport moisture away from the skin, which leaves the user dryer and more comfortable and less likely to suffer from skin problems and ulceration. The soft, lightweight material conforms perfectly to the user's body shape providing a comfortable fit.

For more information visit www.safehip.com

RECOMMENDED SIZES

SIZE	CM	IN
XS	65-85	26-34
S	75-95	30-38
M	90-110	36-44
L	100-120	40-48
XL	110-140	44-56
2XL	120-150	48-60

Measurement at the hip

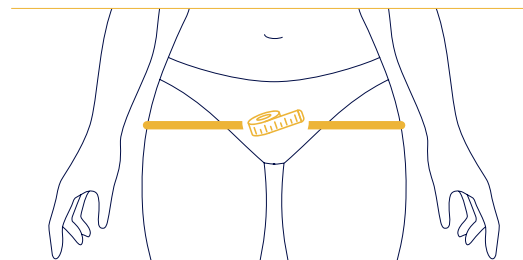

CONTACT INFORMATION

Corsinel Tytex A/S
Marsvej 6
DK-7430 Ikast
mail@tytex.com

COMPOSITION

58% Cotton
36% Polyamide
6% Elastane

AVAILABLE IN SIZES



CARE HANDLING



The textile shields allow easy washing and drying of the SAFEHIP® AirX hip protectors. Best durability achieved at 60° C / 140° F washing and tumble dried at the highest possible drying temperatures without damaging the shields.