

DISCREET AND COMFORTABLE
HIP PROTECTION

SAFEHIP[®] AIR X DISCREET



ABOUT SAFEHIP AIR X DISCREET

> The protective shields in SAFEHIP[®] AirX Discreet are based on the horseshoe-shape, which not only reduces the energy of the impact of the fall, but also diverts it into the surrounding soft tissue.

> The SAFEHIP[®] AirX Discreet hip protectors have a unique and thin design which makes them discreet to wear – thereby providing ultimate discretion as well as protection. The pant and the 100% textile shields are moisture transporting and breathable which makes them very comfortable to wear.

FEATURES

- Comfortable to wear and invisible to others
- 100% textile all the way through
- 100% breathable
- Considerably more comfortable than conventional hip protectors with rigid non-textile shields
- Easy to take on and off
- Easy to wash and dry
- Designed to resist multiple falls



HIP FRACTURE PREVENTION



SEAMLESS KNITTING



BIOMECHANICALLY TESTED



FACT-BASED SIZING

EUROPEAN-BASED PRODUCTION
FOR OPTIMAL QUALITY,
SUPPLY CHAIN EFFICIENCY
AND FLEXIBILITY



STANDARD
100
6234-7001
DTI

PRODUCT RANGE

SAFEHIP[®] AIR X DISCREET



SAFEHIP[®] AIRX DISCREET
FEMALE WHITE

Art. no. 3604

SAFEHIP[®] AIRX DISCREET
FEMALE BLACK

Art. no. 3604

SAFEHIP[®] AIRX DISCREET
MALE WHITE

Art. no. 3612

SAFEHIP[®] AIRX DISCREET
MALE BLACK

Art. no. 3612

PRODUCT DESCRIPTION

The designs of SAFEHIP[®] AirX Discreet hip protector pants ensure optimum protection and come in two models – a male and a female. The pants are flexible and comfortable to wear, not too warm, and may be used with small and light incontinence pads as fixation pants.

Both versions have been designed to perform like regular underwear and can be worn over normal underwear as well.

The SAFEHIP[®] AirX Discreet pants are also very durable and can be washed at 60°C/140°F and tumble dried.

WHAT ARE THE BENEFITS OF AIR X[®] SHIELDS?


AirX is a three dimensional spacer fabric which is ideally suited to absorb impacts and cushion the body from falls. One of the unique features of the material is that it has the ability to effectively transports moisture away from the skin, helping to keep the user cooler than any other hip protectors on the market, which leaves the user dryer and more comfortable and less likely to suffer from any skin problems and ulceration. The soft, lightweight material follows the contours of the body, thus the flexible AirX conforms perfectly to the user's body shape - resulting in improved user comfort.

For more information visit www.safehip.com

RECOMMENDED SIZES

SIZE	CM	IN
S	75-95	30-38
M	90-105	36-42
L	100-115	40-46
XL	110-130	44-52
2XL	125-150	50-60

Measurement at the hip

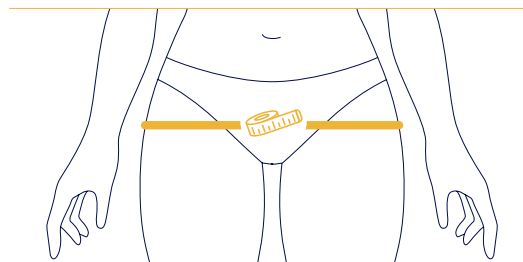

CONTACT INFORMATION

Corsinel Tytex A/S
Marsvej 6
DK-7430 Ikast
mail@tytex.com

COMPOSITION

Shield
89% Polyester
11% Elastane

Seamless briefs
94% Polyamide
6% Elastane

AVAILABLE IN SIZES



CARE HANDLING



The textile shields allow easy washing and drying of the SAFEHIP[®] AirX Discreet hip protectors. They have high durability and can be washed at temperatures up to 60°C and tumble dried without damaging the shield.